

Introducing

Topi-CLICK[®] MICRO[™]

Revolutionary Way to Accurately & Easily Micro-Dose
Topical Medications for Humans & Pets



**Targeted Topical
Micro-Dosing**

Easy Screw Cap

Choice of three
applicators

Standard Luer
Lock Connection

REFILL Indicator

Innovative
TRIAD Shape
Allows Precision
Application
Control

Three Visual
CLICK Points

Each turn/CLICK of
Base dispenses 0.05 mL
of topical

1 CLICK = 0.05 mL
2 CLICKs = 0.10 mL
3 CLICKs = 0.15 mL



Base

The Topi-CLICK[®] Micro[™] Difference

Accurate & Consistent Dose[™]

- Topi-CLICK[®] dispensers proven most accurate among others tested by *ARL
- Eliminates user-variability
- Consistent automated dosing throughout treatment

Automated-Metered Dosing

- Audible, Tactile and Visual CLICKs indicate measure of dose at each point of the unique triad shaped base
- One CLICK = 0.05 mL

Targeted Topical Therapy

- Innovative triad shaped dispenser allows for precision handling during application
- Applicator reduces the risk of hand transference to others
- Offers applicator selections for customized dosing

Select one of 3
applicators



Soft
Angled
Applicator[™]



Pin Point
Applicator[™]



Rounded
Applicator[™]

Viscosity Assurance[™]

- Accuracy NOT affected by varying formulations
- Accuracy NOT affected by varying bases

Dispense/CLICK: 0.05 mL

Delivery Volume: 1 mL to 9.75 mL

Approx. Residual: 0.15 mL

Made of FDA grade materials | Phthalate free | UV block and is Recyclable ♻️.

Topi-CLICK[®] MICRO[™]


User Instructions



- ✓ Automated-metered topical dosing provides an accurate and consistent dose every time
- ✓ Audible, Tactile and Visual CLICKs indicate measure of dose (no need to manually mark dose)
- ✓ Innovative TRIAD shape allows for precision handling
- ✓ Built-in REFILL Badge helps patients reorder medication on-time



STEP 1:
Unscrew the cap of the Topi-CLICK[®] Micro[™] and set aside for recapping after use.

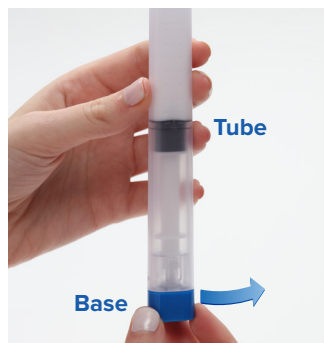

STEP 2:
Holding the Micro[™] in an upright position, turn/CLICK the **Base** the number of CLICKs or mLs specified on the prescription label. Use the Audible, Tactile and Visual CLICK points of the triad shaped **Base** as dosing indicators. Be sure to align **Base** with the **Triad Tube** to ensure dose is fully dispensed.



STEP 3:
Apply topical medication with built-in applicator as directed by physician or pharmacist. **Continue rubbing until preparation is completely absorbed into the skin.**
IMPORTANT: The built-in applicator lessens the risk of transference to unintended areas. No need to use your fingertips to apply medication.



STEP 4:
Screw-on cap after each dose to prevent evaporation. Be sure cap is secured completely.



The Refill Badge[™] reminder system helps patients reorder in a timely manner.

When the **TOP** black line of the **Plunger** reaches the embossed line below the word "REFILL" approximately 32 CLICKs are left. Approximately 12 CLICKs remain at the line above the word "REFILL".

Note: Follow instructions on prescription label before use. Topi-CLICK[®] Micro[™] is not re-fillable.

Made of FDA grade materials | Phthalate free | UV block and is Recyclable ♻️.

Topi-CLICK[®] MICRO[™] Filling Instructions

Micro Quick-Fill Technology[™]
Makes Filling Fast & Simple


1. BEFORE FILLING

Prime syringe or EMP Jar to insure only medication (no air) is loaded into the Micro[™].

Graduated Lines & Top Line of Plunger Indicate Fill Volume.

Design prevents air from entering tip when filling.

Standard Luer Lock to "female". No need for extra connector.

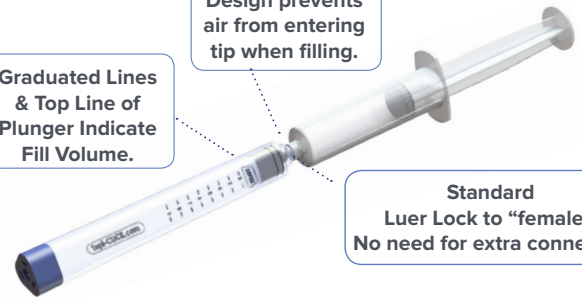

2. FILLING

Note: Graduated Lines are for filling only and not for patients. Top black line of the plunger indicates Fill Volume. The embossed Refill Badge helps patients reorder medication in a timely manner.

Syringe Tip-to-Tip Filling:
Screw Luer syringe to the Luer opening of the Micro[™]. Requires no extra Luer Lock adaptor.



Tip-to-tip Filling

EMP Jar Quick-Fill Adaptor[™] Filling:
Screw-on EMP Jar Quick-fill Adaptor[™]. Fill topical medications directly from mixing jar with this added Micro[™] accessory.



EMP Jar Quick-Fill Adaptor[™]

3. APPLICATOR TIP ATTACHMENT

Select Applicator Tip and press on by pushing at a 90 degree angle (important) against a hard surface.



Select one of 3 applicators



Soft Angled Applicator[™]



Pin Point Applicator[™]



Rounded Applicator[™]



REFILL BADGE

Note: Soft Angled Applicator[™] placement requires tip to be aligned where lower slant of tip is on the same side as Refill Badge of the Micro[™] tube.

4. PRIMING MICRO[™]

Applicator Tip must be attached before priming. Prime the Micro[™] by turning/CLICKing the base until contents are dispensed, then wipe off. Be sure to align Micro[™] base with the triad tube. Approximate Priming: 0.10 mL or 2 CLICKS.



Base

5. CAP PLACEMENT

Screw-on Cap securely to prevent evaporation.



6. LABEL APPLICATION

Apply label to the bottom side of the Topi-CLICK[®] Micro[™], opposite the "REFILL Badge" indicator.



Dispense/CLICK: 0.05 mL
Delivery Volume: 1 mL to 9.75 mL
Approx. Residual: 0.15 mL



DoseLogix Laboratory Test Results Prototype 170706-1

Product performed per design specifications. Device is both highly accurate and precise, demonstrating dosage tight variance as well as an average dose being +1.29% greater than the rated. Product performance variance was reduced for two CLICKS versus one CLICK.

Product prototype 170706-1 qualifies inspection and is passed.

Accurate Micro Dosing Designed for Patients and Pets

Topi-CLICK® Micro Lifetime Dosing Accuracy (per 2 CLICKS)

Average Dose	0.101 mL
Standard Deviation	0.004 mL
Coefficient of Variation	4.46%
Average Dose to Rated Volume (2 CLICKS or 0.10 mL)	+1.29%

Dispense/CLICK: 0.05 mL
Dispenser Volume: Up to 9.75 mL
Delivery Volume: Up to 9 mL

Approx. Residual: 0.15 mL
Approx. Priming: 0.10 mL or 2 CLICKS

*Volumes are approximate. As required for compounded prescription pharmaceuticals dispensed in volume dispensers; verify the weight to volume conversions thoroughly.

Made of FDA grade materials | Phthalate free | UV block and is Recyclable 



DoseLogix Laboratory Test Results Prototype 170706-1

Product performed per design specifications. Device is both highly accurate and precise, demonstrating dosage tight variance as well as an average dose being +1.29% greater than the rated. Product performance variance was reduced for two CLICKS versus one CLICK.

Product prototype 170706-1 qualifies inspection and is passed.

Accurate Micro Dosing Designed for Patients and Pets

Topi-CLICK® Micro Lifetime Dosing Accuracy

Average Dose	0.051 mL
Standard Deviation	0.004 mL
Coefficient of Variation	8.1%
Average Dose to Rated Volume (0.050 mL)	+1.29%

Dispense/CLICK: 0.05 mL
Dispenser Volume: Up to 9.75 mL
Delivery Volume: Up to 9 mL

Approx. Residual: 0.15 mL
Approx. Priming: 0.10 mL or 2 CLICKS

*Volumes are approximate. As required for compounded prescription pharmaceuticals dispensed in volume dispensers; verify the weight to volume conversions thoroughly.

Made of FDA grade materials | Phthalate free | UV block and is Recyclable