

Topi-CLICK[®]

Topi-CLICK[®] MICRO[®]

Revolutionary way to accurately & easily micro-dose topical medications for humans & pets



0.05 mL
Per
CLICK



NEW
Colors Available
Topi-CLICK Micro[®] 9 mL

Easy Screw Cap



Choice of three
applicators




Standard Luer
Lock Connection



REFILL Indicator



Innovative
TRIAD Shape
Allows Precision
Application
Control

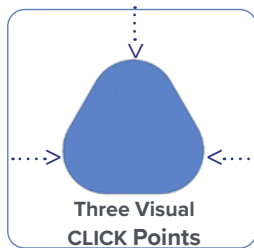

Each turn/CLICK of
Base dispenses 0.05 mL
of topical

1 CLICK = 0.05 mL
2 CLICKS = 0.10 mL
3 CLICKS = 0.15 mL



Base

The Topi-CLICK Micro[®] Difference

Accurate & Consistent Dose[™]

- Topi-CLICK Micro[®] is the most accurate among those in a patient simulation test.
- Eliminates user-variability
- Consistent automated dosing throughout treatment

Automated Metered Dosing

- Audible, Tactile and Visual CLICKs indicate measure of dose at each point of the unique triad-shaped base (One CLICK = 0.05 mL)

Targeted Topical Therapy

- Innovative triad-shaped dispenser allows for precision handling during application
- Applicator reduces the risk of hand transference to others
- Offers applicator selections for customized dosing

Select one of 3
applicators



Universal
Rounded
Applicator[™]



Soft
Angled
Applicator[™]



Pin Point
Applicator[™]

Viscosity Assurance[™]

- Accuracy NOT affected by varying formulations
- Accuracy NOT affected by varying bases

Topi-CLICK Micro[®] 9 mL

Dispense/CLICK: 0.05 mL

Dispenser Volume: Up to 9.75 mL

Total Delivery Volume: Up to 9.45 mL

Approx. Residual: 0.15 mL

Approx. Priming: 0.15 or 3 CLICKs

Made of FDA grade materials | Phthalate free | blocks UV and is Recyclable ♻️.

DOSELOGIX

Copyright © 2020, DoseLogix. All rights reserved.

DoseLogix.com | 877.870.8448 | info@DoseLogix.com

Request FREE samples for patient
trials at DoseLogix.com/samples/

Topi-CLICK[®]

Topi-CLICK[®] MICRO[®] User Instructions




Topical Micro-dosing Application

- ✓ Automated, metered topical dosing provides an accurate and consistent dose every time
- ✓ Audible, Tactile and Visual CLICKs indicate measure of dose (no need to manually mark dose)
- ✓ Innovative triad-shape allows for precision handling
- ✓ Built-in REFILL Badge helps patients reorder medication on-time



STEP 1:
Unscrew the cap of the Topi-CLICK Micro[®] and set aside for recapping after use.

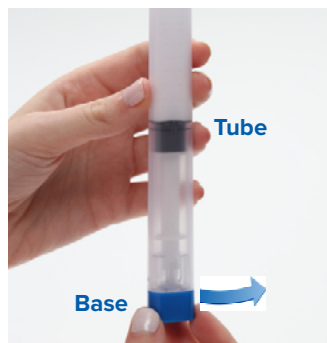

STEP 2:
Holding the Micro[®] in an upright position, turn/CLICK the **Base** the number of CLICKs or mLs specified on the prescription label. Use the Audible, Tactile and Visual CLICK points of the triad-shaped **Base** as dosing indicators. Be sure to align **Base** with the **Triad Tube** to ensure dose is fully dispensed.



STEP 3:
Apply topical medication with built-in applicator as directed by physician or pharmacist. **Continue rubbing until preparation is completely absorbed into the skin.**
IMPORTANT: The built-in applicator lessens the risk of transference to unintended areas. No need to use your fingertips to apply medication.



STEP 4:
Screw on cap after each dose to prevent evaporation. Be sure cap is secured completely.

Innovative TRIAD Shape Allows Precision Application Control

Each turn/Click of the base dispenses 0.05 mL of medication accurately.

REFILL Badge™



REFILL Line

Plunger

The Refill Badge™ reminder system helps patients reorder in a timely manner.

When the **TOP** black line of the **Plunger** reaches the embossed line below the word "REFILL" approximately 32 CLICKs are left. Approximately 12 CLICKs remain at the line above the word "REFILL".

Note: Graduated Marks are for pharmacy filling purposes only. Follow instructions on prescription label before use. Topi-CLICK Micro[®] is not re-fillable.

Made of FDA grade materials | Phthalate free | blocks UV and is Recyclable ♻️.

DOSELOGIX

Copyright © 2020, DoseLogix. All rights reserved.

DoseLogix.com | 877.870.8448 | info@DoseLogix.com

Request FREE samples for patient trials at DoseLogix.com/samples/



Topi-CLICK

MICRO[®]

DoseLogix Laboratory Test Results

Product performed per design specifications. Device is both highly accurate and precise, demonstrating a tight variance and an average dose being 98.7% of the rated dose. Product performance variance was reduced for two CLICKs versus one CLICK.

Qualifies inspection and is passed.

Accurate Micro-dosing Designed for Patients and Pets

Topi-CLICK Micro[®] Lifetime Dosing Accuracy

(1 CLICK or 0.05 mL)

(2 CLICKS or 0.10 mL)

Average Dose	0.049 mL	0.099 mL
Standard Deviation	0.003 mL	0.004 mL
Coefficient of Variation	5.7%	3.8%
Average Dose to Rated Volume	98.7%	98.7%

Dispense/CLICK: 0.05 mL

Dispenser Volume: Up to 9.75 mL

Total Delivery Volume: Up to 9.45 mL

Approx. Residual: 0.15 mL

Approx. Priming: 0.15 mL or 3 CLICKS

*Volumes are approximate. As required for compounded prescription pharmaceuticals dispensed in volume dispensers; verify the weight to volume conversions thoroughly.

Made of FDA grade materials | Phthalate free | blocks UV and is Recyclable 



Copyright © 2020, DoseLogix. All rights reserved.

DoseLogix.com | 877.870.8448 | info@DoseLogix.com

Request FREE samples for patient
trials at DoseLogix.com/samples/