



- > Appeals to Safety, Not to Children
- > Push-and-Turn Security Cap
- > DoseLock™ Base with Integrated Safety Tabs
- > Senior-Friendly Design

NEW!

The Topi-CLICK® 35 Child Resistant topical dosing dispenser delivers the same trusted performance and dosing accuracy¹ as the standard Topi-CLICK® 35, with the added protection of a certified-child resistant design for enhanced safety.

The dosing applicator can be used to dispense a total volume of up to 35 mL of a topical preparation, depending on the volume-to-weight ratio of the formulation. In addition to providing accurate and repeatable volumetric dosing, the integrated applicator design helps ensure that topical medications are applied only to the intended area, thereby reducing the risk of improper dosing and contamination.

Certified Child Resistance

*Meets U.S. Consumer Product Safety Commission (CPSC) criteria for child-resistance under 16 CFR §1700.20.

Push-and-Turn Security Cap

Deliberately familiar locking cap with push-down-and-turn design is certified to resist access by children while remaining intuitive for adult users.

Topi-CLICK DoseLock™ Base with Integrated Safety Tabs

Pinch-grip dispensing mechanism requires coordinated adult-level dexterity for activation, preventing unintended operation by children.

Two-Stage Access Protection

Medication is inaccessible unless both locking cap is removed and pinch-grip base is actively engaged—providing layered safety.

Built-In User Guidance

Clear instructional text and directional arrows printed on the cap and base to ensure easy and reliable operation by adult users.

Appeals to Safety, Not to Children

Universal gray cap and base reduce visual appeal to children, while multiple cap color options still allow for subtle distinction between medications.

* Certified Child Resistance. Meets U.S. Consumer Product Safety Commission (CPSC) criteria for child-resistance under 16 CFR §1700.20. (only when fully assembled as specified).

¹ Request independent accuracy study at https://topi-click.com/accuracy_study/.



- Accuracy assures patients request refills on-time, every time.
- Built-in applicator for hands-free application, which lessens the risk of cross-contamination.
- Multiple applicator color options for differentiating medications.
- Audible, visual and tactile CLICKs make it easy for patients to accurately dose medications.
- Has a Refill Badge™ for prompting patients to act for on-time refills.



Child-Resistant Cap

.25 mL
Per CLICK

DoseLock™
Safety Tab

Applicator Pad and
Cap is available in
Blue, Gray, Pink or White



Topi-CLICK | TEAM Technologies
300 Parkbrooke Place, Suite 210, Woodstock, GA 30189
877-870-8448 | info.TopiCLICK@teamtech.com
www.Topi-CLICK.com



Pharmaceutical
Delivery
Devices

INSTRUCTIONS FOR USING TOPI-CLICK[®] 35 CHILD RESISTANT

STEP 1: Unlock and Remove Child Resistant Cap



Read instructions on prescription label.

Remove the *Child-Resistant Cap. Place palm over the cap and press downward firmly while turning counterclockwise. Once unlocked, lift the cap off the dispenser.

Tip: The top of the Child-Resistant Cap includes the instructional text "PUSH DOWN & TURN" as well as directional arrows labeled "OPEN" and "CLOSE" as a convenient visual reminder of the proper technique for cap opening and closure.

STEP 4: Apply Medication Using Built-in Applicator



Apply topical medication by rubbing the area of skin as directed by your pharmacist using the **built-in applicator**. Continue rubbing until topical is completely absorbed into the skin.

IMPORTANT: The Built-in Applicator lessens the risk of cross-contamination and ensures the full dose is applied to the intended area.

STEP 2: Locate Two Safety Tabs-each labeled PUSH on DoseLock[™] Base



Locate the two safety tabs - each labeled "PUSH"- on the base of dispenser.

Simultaneously press two DoseLock[™] Safety Tabs (i.e., tabs labeled "PUSH") of base using grip with one hand while holding dispenser tube firmly in other hand.

STEP 5: Reinstall and Lock Child-Resistant Cap



Replace cap after each dose to prevent evaporation. Align the Child-Resistant Cap over the applicator and press down while turning clockwise until the cap is securely locked.

STEP 3: Turn DoseLock[™] Base to Dispense Medication



While maintaining pinch grip, rotate the DoseLock[™] base to deliver the number of CLICKs specified in prescription.

STEP 6: Monitor the Refill Indicator



The bottom edge of the plunger is the indicating line; when it reaches the embossed line below the word "Refill" approximately 32 CLICKs are left. Approximately 12 CLICKs remain at the line above the word "Refill".

*Certified Child Resistance. Meets U.S. Consumer Product Safety Commission (CPSC) criteria for child-resistance under 16 CFR §1700.20. (only when fully assembled as specified.)

Topi-CLICK[®] 35 is the most accurate topical dosing dispenser among those compared in an independent study titled Topical Metered Dosing Dispenser Performance Evaluation by Analytical Research Laboratories. Request study at https://topi-click.com/accuracy_study/.

PRECAUTIONS

- Store in a dry place at room temperature.
- Avoid exposure to sunlight and other bright lights.
- When traveling in an airplane or to a different elevation, seal the applicator in a zip lock bag to prevent soiling from any leakage due to pressure changes.



Topi-CLICK | TEAM Technologies
 300 Parkbrooke Place, Suite 210, Woodstock, GA 30189
 877-870-8448 | info.TopiCLICK@teamtech.com
www.Topi-CLICK.com



**Pharmaceutical
 Delivery
 Devices**

FILLING THE APPLICATOR



STEP 1: Fill the Dispenser Tube

Fill the Topi-CLICK[®] 35 *Child Resistant dispenser with the desired amount of topical formulations using a Topi-CLICK[®] RapidFill™ System, spatula, EMP jar, or nozzle.



STEP 2: Settle the Formula (if needed)

If needed while filling, gently tap the base of the dispenser on your palm, counter, or a Tap-n-CLICK™ device to help settle the topical formula and remove air pockets.




STEP 3: Install the Child-Resistant Applicator

Align the internal key tabs of the applicator with the receiving slots at the top of the dispenser tube. Press down until the applicator snaps into place.



Locking Feature

NOTE: Only use the Child-Resistant Applicators that are supplied with the Topi-CLICK[®] 35 Child Resistant dispenser product. Always verify that the appropriate, child-resistant applicator is installed by confirming presence of the distinct, raised locking features located around the outer edge of the applicator.

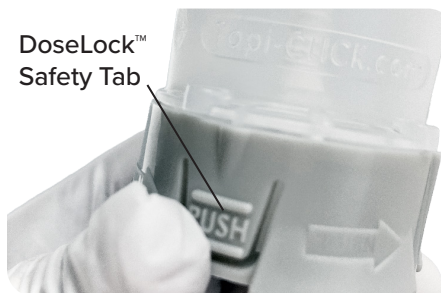

STEP 4: Prime the Dispenser

Simultaneously press two DoseLock™ Safety Tabs (i.e., tabs labeled "PUSH") of base using pinch grip with one hand while holding dispenser tube firmly in other hand. While maintaining this pinch grip, slowly rotate the base until a small amount of medication is dispensed from the applicator.



Wipe excess topical from the applicator after priming is complete.

NOTE: Priming typically requires approximately 2 CLICKs (0.50 mL).



STEP 5: Install the Child-Resistant Locking Cap

Place the outer cap on top of the dispenser. Push down firmly while turning clockwise to engage and lock the child-resistant mechanism.



NOTE: Verify the cap is securely locked by confirming it cannot be removed without using the push-down-and-turn action.

**Certified Child Resistance. Meets U.S. Consumer Product Safety Commission (CPSC) criteria for child-resistance under 16 CFR §1700.20. (only when fully assembled as specified.)*



STEP 6: Apply the Label

Apply a prescription label to the tube.

NOTE: Do not cover the Refill badge.



Topi-CLICK | TEAM Technologies
 300 Parkbrooke Place, Suite 210, Woodstock, GA 30189
 877-870-8448 | info.TopiCLICK@teamtech.com
www.Topi-CLICK.com



**Pharmaceutical
 Delivery
 Devices**



Topi-CLICK® 35 Child Resistant: 35 mL, 0.25 mL per CLICK

DESCRIPTION: The Topi-CLICK® 35 Child Resistant topical dosing dispenser delivers the same trusted performance and dosing accuracy as Topi-CLICK® 35, with the added protection of a certified-child resistant design for enhanced safety.

Topi-CLICK® 35 is designed to overcome the various limitations associated with pumps, syringes, tubes, tubs and jars and can be used to dispense approximately 5 mL to 35 mL of a topical preparation depending on the volume to weight ratio of base and formula used.

Topi-CLICK® makes it easy to accurately measure the dose and apply topical medications to the intended area, thereby decreasing the risk of improper dosing and contamination while improving compliance and outcomes of therapy.

FEATURES: Topi-CLICK® 35 Child Resistant

- Child Resistance - Meets U.S. Consumer Product Safety Commission (CPSC) criteria for child-resistance under FDA 16 CFR §1700.20 (only when fully assembled as specified.)
- Each CLICK or turn of the base moves the plunger approximately 0.25 mL of volume space, thereby dispensing 0.25 mL of cream on average.
- The most common number of CLICKs used for one dose is 2 CLICKs or 0.50 mL of volume space.
- Audible, visual and tactile CLICK registrations for user recognition.
- Lower cost with simplified labeling – labels as fast and easily as a prescription bottle.
- Built-in applicator to prevent medication transference to unintended areas.
- “Refill” reminder badge helps patients to reorder in a timely manner.
- Designed to resist accidental discharge in purse, pocket or luggage.
- Multiple applicator color options for differentiating medications.

TECHNICAL SPECIFICATIONS: Topi-CLICK® 35 Child Resistant

Dispenser Volume:	Up to 37 mL
Delivery Volume:	Up to 34 mL
Dispensed per CLICK:	0.25 mL
Approximate Residual:	2.5 mL
Approximate Priming:	0.50 mL or 2 CLICKs
Doses Provided:	When filling with 35 mL, the Topi-CLICK® will provide approximately 32 mL or 64 doses (at 2 CLICKs per dose or 128 CLICKs total).

- Registered FDA device, Class 1 Dispenser.
- FDA Registration # is 3007142798.
- Materials compliant with FDA 21.177 code, which is the requirement for United States Pharmacopeia (USP <661>).
- Made of FDA food grade materials | BPA and Phthalate free.
- Topi-CLICK® 35 Child Resistant is USP Compliant to USP671 (250nm to 450 nm) and blocks UVA, UVB and Violet Light.
- Recycling Code is #5 ♻️. Visit <http://Topi-CLICK.com/recycling-your-topi-click> for more information on recycling Topi-CLICK®.

NOTE: Volumes are approximate. As required for compounded prescription pharmaceuticals dispensed in volume dispensers, verify the weight to volume conversions thoroughly. See dispensing instructions at <http://Topi-CLICK.com/downloads-how-it-works>.

One time use only.

Topi-CLICK® 35 is the most accurate topical dosing dispenser among those compared in an independent study titled *Topical Metered Dosing Dispenser Performance Evaluation by Analytical Research Laboratories*. Request study at https://topi-click.com/accuracy_study/.



Topi-CLICK | TEAM Technologies
300 Parkbrooke Place, Suite 210, Woodstock, GA 30189
877-870-8448 | info.TopiCLICK@teamtech.com
www.Topi-CLICK.com



**Pharmaceutical
Delivery
Devices**