




- ✓ Flexible Easy-press Buttons™ to eject full troches.
- ✓ Dispense directly to patient with disposable mold design.
- ✓ Easy-to-eject and reinsert for quality control checks.
- ✓ Heat resistant tray means less warping when filling or removing excess base with heat gun.
- ✓ Large label surface on cover.
- ✓ No need for extra sleeve covers.

Filling Instructions



STEP 1:
Place the Troche-CLICK™ tray on a clean, level surface. Fill the Troche-CLICK™ with mixture.

Filling techniques may include: pouring directly from a glass beaker, transferring mixture using a syringe or transfer pipette.

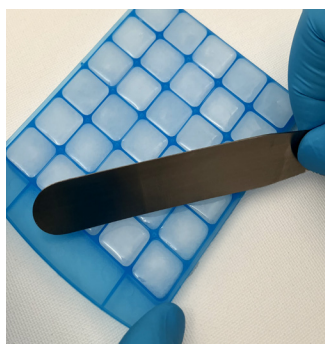


STEP 2:
Remove excess troche material by scraping over the top of the troche lozenges with a clean spatula.

The PEG base shrinks during cooling. Add additional material to the individual cavities, if required.



STEP 3:
After the troche material has set, a heat gun can be used to remelt the top of the troches to provide a polished appearance to the troches.

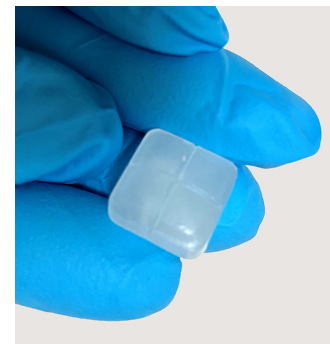

STEP 4:
Troches lozenges can be removed and then replaced to conduct quality checks.



STEP 5:
Slide the troche tray into the cover. Label can now be applied to the cover and lip of the tray.

TECHNICAL SPECIFICATIONS:

Topi-CLICK® Troche-CLICK™

Dispenser Volume: Up to 30 individual 1.0g (PEG 1450) troche lozenges.

Delivery Volume: 1.0g (PEG 1450) per troche (Each troche can be divided into halves or quarters).

Doses Provided: Device can hold thirty 1.0g (PEG 1450) doses.



Topi-CLICK® Troche-CLICK™ Specifications

Topi-CLICK® Troche-CLICK™ : Each cavity calibrated to 1.0 gram with PEG 1450. Scored halves and quarters allow customizable dosing.

DESCRIPTION: Topi-CLICK® Troche-CLICK™ is an innovative way to provide an accurate dose to patients in a troche lozenge. The easy-to-fill Troche-CLICK™ increases pharmacy efficiency and improves the patient experience with medications dispensed directly from the disposable mold design. Troche-CLICK® is available in 4 vibrant colors which allows for easy medication differentiation and quality control.

FEATURES:

- Each cavity will contain a 1.0g (PEG 1450) smooth troche lozenge designed for comfortable usage.
- Each cavity provides score lines for easy division of troches into halves or quarters for greater customizability of dosing.
- Each cavity has a textured, easy eject button for quick removal of troches from the tray.
- Large labeling surface on cover.
- Secondary labeling surface on tray.
- Translucent tray allows for easy view of doses remaining.
- CLICK locking mechanism designed to resist accidental opening.
- Manufactured from heat resistant materials.

TECHNICAL SPECIFICATIONS: Topi-CLICK® Troche-CLICK™

Dispenser Volume: Up to 30 individual 1.0g (PEG 1450) troche lozenges.

Delivery Volume: 1.0g (PEG 1450) per troche (Each troche can be divided into halves or quarters).

Doses Provided: Device can hold thirty 1.0g (PEG 1450) doses.


- Class 1 Medical Device
- Topi-CLICK's FDA Registration # is 3007142798.
- Materials compliant with FDA CFR 21.177 code, which is the requirement for (<USP 661>).
- Made of FDA grade materials | BPA and Phthalate free | Blocks UV.
- Recycling Code is #5. ♻️ Visit <http://https://topi-click.com/recycling-info> for more information on recycling Topi-CLICK®.

***Volumes are approximate. As required for compounded prescription pharmaceuticals dispensed in volume dispensers, verify the weight to volume conversions thoroughly. One time use only.**



- ✓ Flexible Easy-press Buttons™ to eject full troches.
- ✓ CLICK locking mechanism eliminates accidental opening.
- ✓ Translucent tray allows for easy view of doses remaining.
- ✓ Available in 4 colors to aid in medication differentiation.

Flexible Easy-press Buttons™
make it easy to eject full troches.

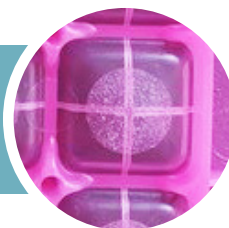

Patient Use Instructions



STEP 1:

Remove troche tray from protective cover.



STEP 2:

Use Easy-press Buttons™ to eject troche lozenge. Dispense troche lozenge dosage as directed by your physician or pharmacist.



STEP 3:

Slide protective cover over tray to activate the CLICK locking mechanism.

Topi-CLICK® Troche-CLICK™ is Available in 4 Colors



Multiple Color Options make it easy for patients to distinguish between their medications.